



Encountering a Sunny day

- 2 Handpieces
- 3 Safety Securing Systems
- 4 Advanced Techonolgies



Category	Sunny	
Energy Type	Radio Frequency	
HP Name	Pulsar	Sunny
Level	1 – 10	0.5 – 8.0
Pulse	1 – 10	6
Mode	X	Soft & Hard
Polar	Mono & Bi	Mono
Cooling	X	ACA
Frequency	1 MHz, 2 MHz	6.78MHz
Tip Type	Quasar Tip, Blazar Tip	Sunny Face, Sunny Eye
RF-Time	100-800 ms	Fix
Impedance Check	Real Time Impedance Check	
Switch	Foot, Hand	
Display	12.1 inch Wide LCD Touch Screen	
Consumption	700 VA	
Dimension	(W)681 mm x (L)467 mm x (H)1415 mm	
Weight	46 (kg)	



Sunny Day, Shining Me

Delivering high-frequency energy into our skin
With SUNNY's unique innovative technologies

Sunny



2 Handpieces

PULSAR

Layer to layer,
Building up energy from
Deep to Upper layers

SUNNY

Sometimes Hard! Sometimes Soft!
Delivering high-density
high-frequency energy

3 Safety Securing Systems

Maximizing Cooling Efficiency Instantaneous Vaporization Technology

Improved cooling efficiency enables
more advantages for pain control management

IVT

Uniform Cooling Effect Advanced Gas Cooling Atomization Technology

Enhancing the safety during treatment

ACA

Real-time Impedance Check

Real-time impedance measurement allows for
consistent energy delivery suited to the skin's condition

RIC

4 Advanced Techonolgies

IMCT

Increase the
contact precision
of 25 patented
high-frequency pins

Individually moving 25 pins
allow precision contact
on various treatment areas



APRT

Sub-pulses up to 10
Interval control system

ADRT

Combination of
Mono-polar & Bi-polar


SWT

Comforting treatment by
Waveform Control system



Before & After

QUASAR 2 treatments

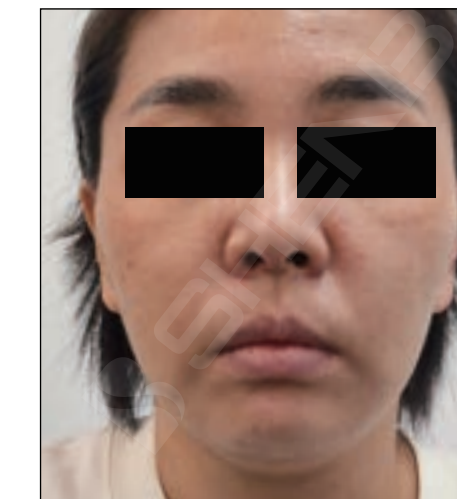

QUASAR 2 treatments



SUNNY FACE 2 treatments



Sunny